



C14H10O4 benzoyl peroxide(BPO)

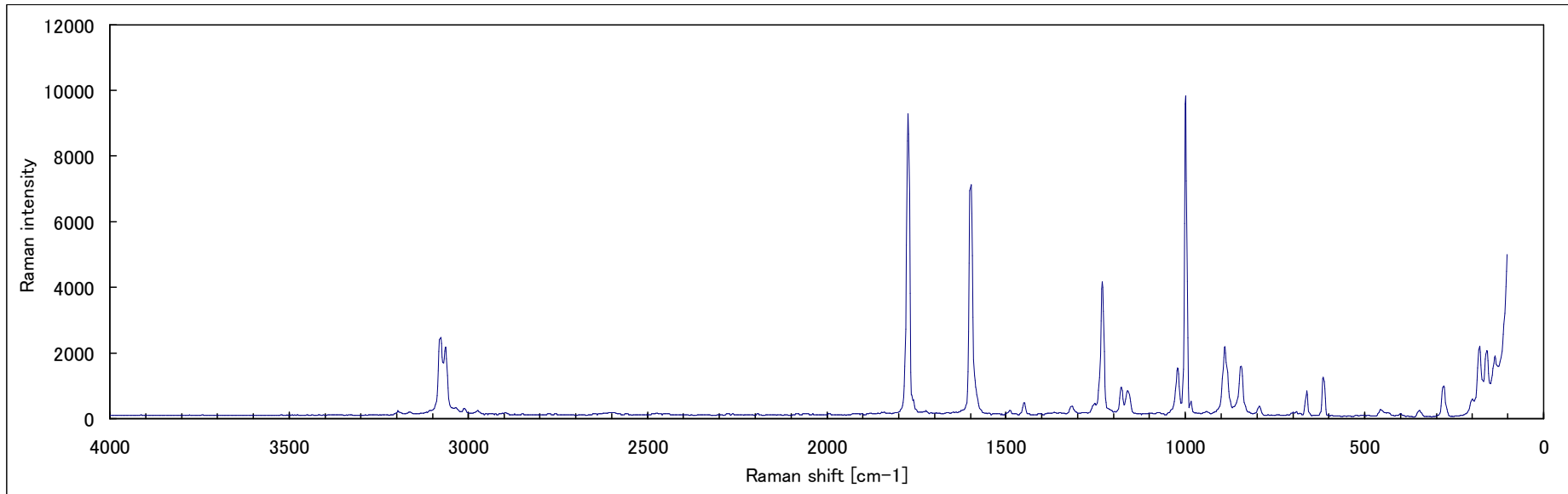

Date Dec 12 2006
 Operator H.F.
 Equipment ALMEGA (Thermo Fisher SCIENTIFIC)
 Laser 532 nm, 25 mW
 Resolution 6.5 – 10.5 cm-1
 Sample solid powder
 Method micro (x50)

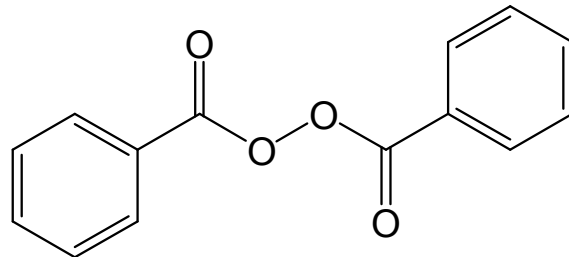

Peak position and Raman intensity

159	1986				
179	2118				
279	900				
615	1162				
661	744				
845	1491				
890	2102				
999	9992				
1021	1451				
1161	743				
1178	873				
1232	4080				
1599	7086				
1774	9207				
3063	2109				
3078	2368				